
Caloric Irrigators Air or Water

The Interacoustics irrigators offer a choice of stimuli, either water or air, for caloric testing during the VNG test protocol.



Air Fx Caloric irrigator

Highlights:

- Cooling and heating ability: 20°C to 50°C
- Lighted and magnified otoscope
- Self-contained cooling ability
- Integrated with the Interacoustics VNG software VN415 / VO425 v. 7.06 or later



Caloric test with the Interacoustics VNG software



Versatile Caloric Air Irrigator

The new AirFx caloric irrigator has many benefits. It offers an otoscopic handle with illumination as well as magnification. In addition, irrigations can be started using the start button on the handle or via the Interacoustics VN415/VO425 software.

The Air Fx is available in 110-130 V and 220-240 V configurations, and is available world-wide.

An integrated water reservoir allows for cool air irrigations down to 20°C regardless of ambient temperature, making it suitable for warm climate use. For warm air irrigations, the Air Fx is capable of heating air up to 50°C.

Temperature Accuracy and Stability is within +/- 0,5°C, and air flow stability is within +/- 0,5 L/min.

Aqua Stim Caloric irrigator



Can be operated from the Interacoustics VNG software VN415 / VO425 v. 7.06 or later



Included external water tank

Highlights:

- 30°C cold and 44°C warm irrigations
- Use with the Interacoustics VNG software or as stand-alone
- LED illumination and start button in handle
- Eliminates the need for drain and faucet in examination room

Versatile Caloric Water Irrigator

The Interacoustics Aqua Stim offers reliable caloric irrigation. The Aqua Stim can be used as a stand-alone device, or you can control it from the Interacoustics VNG software (V. 7.06 or later).

Irrigation is started by simply pressing the start button on the handle one time to initiate the user-defined caloric protocol.

An external water tank has capacity for approximately 14 irrigations, and eliminates the need for a drain in the examination room.

The Aqua Stim is available world-wide and is supplied as a 110-130 V or 220-240 V model.





Interacoustics a/s

Interacoustics is a world leading diagnostic solutions provider in the field of hearing and balance assessment. We help the professional audiologic world reach new milestones through continuous developments and a constant focus on integration and direct customer value.

With you at all times


Whenever you purchase audiological equipment with the Interacoustics brand on it, you are not just guaranteed a quality product but also highly trained support personnel with specialized knowledge. We operate in over 100 countries worldwide through a carefully selected network of distributors and service centres to ensure that you receive total support and service.

Product specifications

All technical and hardware specifications concerning the Air Fx and Aqua Stim Caloric Irrigators can be downloaded from our website.

Read more at
interacoustics.com

Go online to
explore our
full product
range

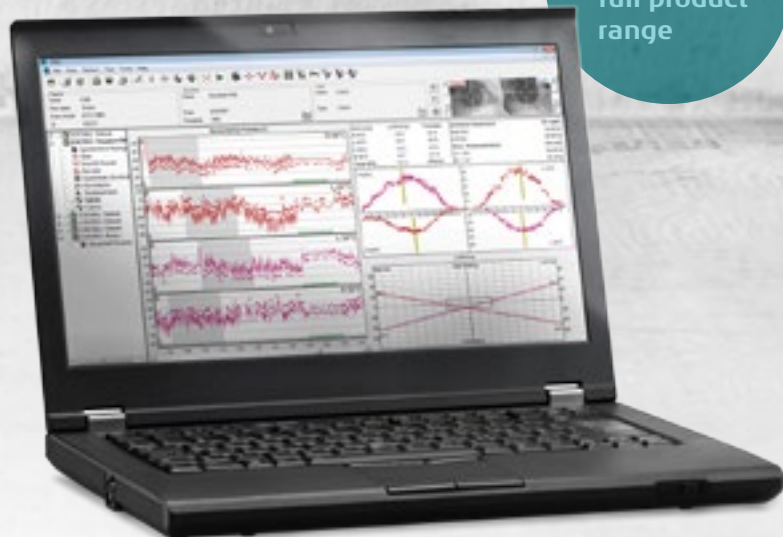

Interacoustics a/s

Drejervaenget 8
DK-5610 Assens

T +45 6371 3555
F +45 6371 3522

info@interacoustics.com
interacoustics.com



Interacoustics